

## Anti-CS1 CAR mRNA-LNP

Ready-to-use lipid nanoparticles

### Order Information

Catalog#	Size	GenPept No.
PM-LNP-0093	200uL	

### Description

CS1 (CD2 subset-1), also known as CD319 (SLAMF7), is a human natural killer (NK) cell receptor encoded by the SLAMF7 gene. This surface antigen is a robust marker of normal and malignant plasma cells in multiple myeloma. Human NK cells express two splice variants of CS1: CS1-S lacks an intracellular domain for activation, while CS1-L contains an intracellular domain and is therefore capable of activating NK cytotoxicity. CS1 is a homophilic receptor, and CS1-CS1 interaction leads to the activation of natural cytotoxicity of NK cells. Both isoforms of CS1 are membrane bound and expressed in NK cells. However, only the CS1-L isoform is expressed in B cells and regulates B cell proliferation and autocrine. This product is designed as a tool for the delivery and expression of human anti-CS1 CAR

### Composition

mRNA-LNPs suspended in PBS (-Ca, -Mg) (pH: 7.0-7.4).

### Translated Protein sequence


Available up request

#### Products and Services

- Monoclonal Antibody
- Bispecific Antibody
- mRNA-LNP
- Hybridoma Sequencing
- CAR-T/NK
- Stable Cell Line

## Anti-CS1 CAR mRNA-LNP

### Ready-to-use lipid nanoparticles

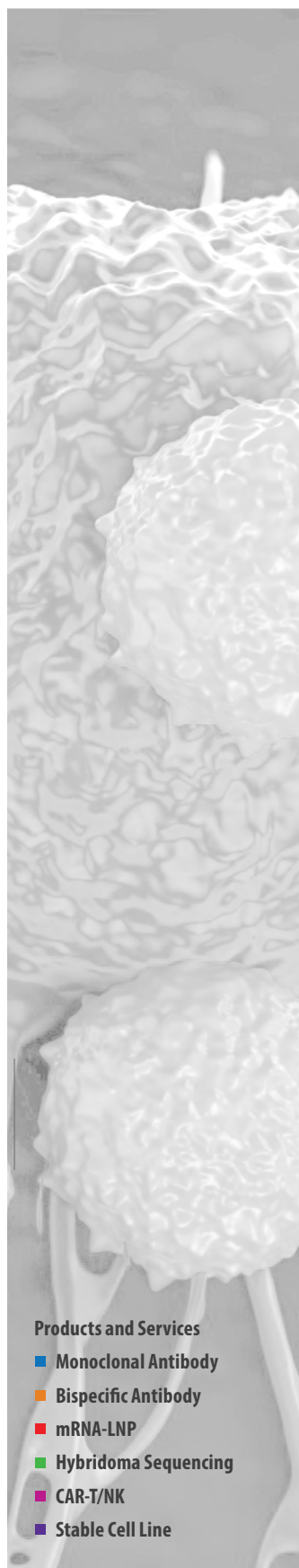

#### Products and Services

- Monoclonal Antibody
- Bispecific Antibody
- mRNA-LNP
- Hybridoma Sequencing
- CAR-T/NK
- Stable Cell Line

#### Storage

Product is delivered on blue ice. Store at 4°C for up to 3 months.

#### Application & Handling

Upon receiving product, briefly pulse spin before opening to ensure product is at bottom of container. It is important not to spin for too long as this may rupture mRNA-LNPs. Do not vortex. Work with mRNA-LNPs on ice and minimize the time that the product spends at room temperature. After handling the product during experiments, return immediately to ice. mRNA-LNP products should only be handled with certified RNase-free reagents and consumables. Use of filtered pipette tips is highly recommended.

#### Safety & Research Disclosure

All ProMab mRNA lipid nanoparticle products are for in vitro research use only. Products are not FDA approved for human use.

#### General Protocol

1. Prior to transfection: Plate 1ml of cells at a density of [1.0E6 cells/ml] in a single well of a 12-well culture plate. Ensure the cells you are using are viable and healthy. Try not to let your cells sit for longer than 5 minutes prior to transfection. Cell clumping at the time of transfection may reduce transfection efficiency. \*Note: If cell clumping occurs, gently pipette your culture up & down to ensure you have a single cell suspension before transfecting.
2. Briefly pipette mRNA-LNP mix up & down to resuspend. Add 20-40ul of the mRNA-LNP product dropwise directly to your 1ml culture. Gently tilt plate back and forth to mix (not necessary if you are using cells which will be immediately placed back into a shaker). Place your transfected cells back into their original culture conditions.
3. Check cell expression by FACS or by using other detection methods at 24hr intervals after transfection. \*Note: This is a generalized protocol for transfection using mammalian suspension culture cells. Transfection volume may be scaled up or down proportionately using the volumes given. HEK-293s cells were grown and transfected in FreeStyle™ F17 Expression Medium (Gibco, Cat#: A1383501), supplemented with GlutaMAX™ (Gibco, Cat#: 35050061), and Poloxamer 188 Non-ionic Surfactant (Gibco, Cat#: 24040032). T-cells were grown and transfected in a culture medium supplemented with 10% FBS (Omega Scientific, Cat#: FB-02). When transfecting cells using mRNA-LNPs, it is typically necessary for the cell culture medium to be supplemented with 10% FBS at the time of transfection. Without this supplement, transfection efficiency will drop significantly. For mRNA-LNP transfection of cells which cannot use FBS in the culture medium, please contact us at (510) 860-4615. Alternative transfection methods are available.

## Anti-CS1 CAR mRNA-LNP

Ready-to-use lipid nanoparticles

### Expression of Human Anti-CS1 CAR in S293 Cells Transfected with PM-LNP-0093

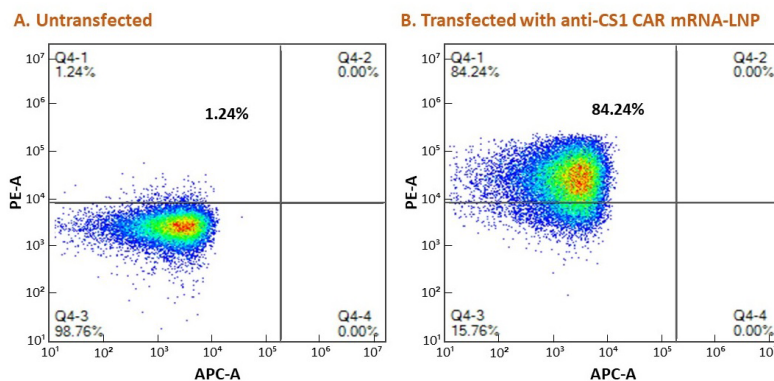


Figure 1. Flow Cytometry. PM-LNP-0093 nanoparticle-treated HEK293S cells express anti-CS1 CAR protein, detected with a labeled anti-F(ab')<sub>2</sub> antibody.

#### Products and Services

- Monoclonal Antibody
- Bispecific Antibody
- mRNA-LNP
- Hybridoma Sequencing
- CAR-T/NK
- Stable Cell Line