

## B7H4-TF mRNA-LNP

Ready-to-use lipid nanoparticles

### Order Information

Catalog#	Size	GenPept No.
PM-LNP-0103	200uL	

### Description

B7H4 (B7 Homolog 4), also known as VTCN1 (V-set domain containing T cell activation inhibitor 1), encoded by the VTCN1 gene, is a protein belonging to the B7 co-stimulatory protein family. These proteins are found on the surface of antigen-presenting cells and interact with ligands bound to receptors on the surface of T cells. B7H4 mRNA is normally expressed in normal tissues, but its protein level is maintained in limited amounts in normal tissues. Studies have shown that high levels of the encoded protein are associated with tumor progression, and in many tumor tissues, B7H4 is abnormally expressed. This suggests that B7H4 expression can be tightly regulated at the post-transcriptional level in normal tissues, but this regulatory mechanism may be suppressed or altered in the tumor environment.

### Composition

mRNA-LNPs suspended in PBS (-Ca, -Mg) (pH: 7.0-7.4).

### Translated Protein sequence

Available up request

#### Products and Services

- Monoclonal Antibody
- Bispecific Antibody
- mRNA-LNP
- Hybridoma Sequencing
- CAR-T/NK
- Stable Cell Line

## B7H4-TF mRNA-LNP

### Ready-to-use lipid nanoparticles

#### Storage

Product is delivered on blue ice. Store at 4°C for up to 3 months.

#### Application & Handling

Upon receiving product, briefly pulse spin before opening to ensure product is at bottom of container. It is important not to spin for too long as this may rupture mRNA-LNPs. Do not vortex. Work with mRNA-LNPs on ice and minimize the time that the product spends at room temperature. After handling the product during experiments, return immediately to ice. mRNA-LNP products should only be handled with certified RNase-free reagents and consumables. Use of filtered pipette tips is highly recommended.

#### Safety & Research Disclosure

All ProMab mRNA lipid nanoparticle products are for in vitro research use only. Products are not FDA approved for human use.

#### General Protocol

1. Prior to transfection: Plate 1ml of cells at a density of [1.0E6 cells/ml] in a single well of a 12-well culture plate. Ensure the cells you are using are viable and healthy. Try not to let your cells sit for longer than 5 minutes prior to transfection. Cell clumping at the time of transfection may reduce transfection efficiency. \*Note: If cell clumping occurs, gently pipette your culture up & down to ensure you have a single cell suspension before transfecting.
2. Briefly pipette mRNA-LNP mix up & down to resuspend. Add 20-40ul of the mRNA-LNP product dropwise directly to your 1ml culture. Gently tilt plate back and forth to mix (not necessary if you are using cells which will be immediately placed back into a shaker). Place your transfected cells back into their original culture conditions.
3. Check cell expression by FACS or by using other detection methods at 24hr intervals after transfection. \*Note: This is a generalized protocol for transfection using mammalian suspension culture cells. Transfection volume may be scaled up or down proportionately using the volumes given. HEK-293s cells were grown and transfected in FreeStyle™ F17 Expression Medium (Gibco, Cat#: A1383501), supplemented with GlutaMAX™ (Gibco, Cat#: 35050061), and Poloxamer 188 Non-ionic Surfactant (Gibco, Cat#: 24040032). T-cells were grown and transfected in a culture medium supplemented with 10% FBS (Omega Scientific, Cat#: FB-02). When transfecting cells using mRNA-LNPs, it is typically necessary for the cell culture medium to be supplemented with 10% FBS at the time of transfection. Without this supplement, transfection efficiency will drop significantly. For mRNA-LNP transfection of cells which cannot use FBS in the culture medium, please contact us at (510) 860-4615. Alternative transfection methods are available.

#### Products and Services

- Monoclonal Antibody
- Bispecific Antibody
- mRNA-LNP
- Hybridoma Sequencing
- CAR-T/NK
- Stable Cell Line

## B7H4-TF mRNA-LNP

Ready-to-use lipid nanoparticles

### Expression of Human B7H4-TF in S293 Cells Transfected with PM-LNP-0103

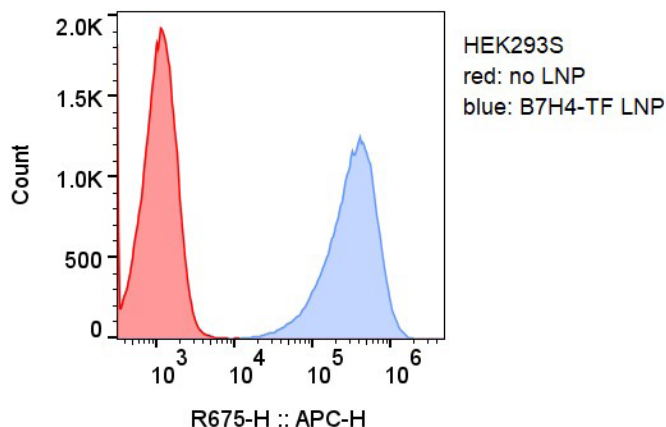


Figure 1. Flow Cytometry. PM-LNP-0103 nanoparticle-treated HEK293S cells express TF-tagged B7H4 protein, detected with a labeled anti-TF antibody.

### Secretion of Human B7H4-TF by S293 Cells Transfected with PM-LNP-0103

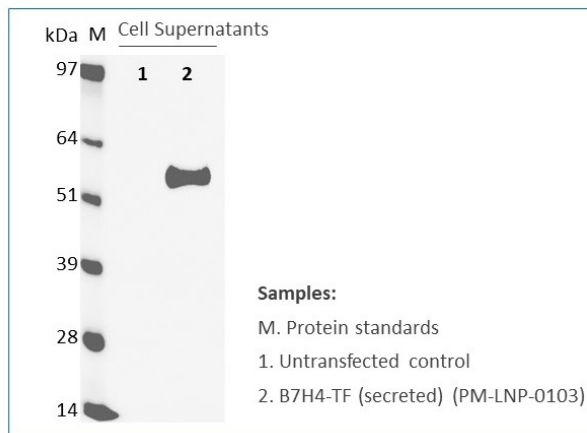


Figure 2. Western Blot. Medium collected from PM-LNP-0103 nanoparticle-treated HEK293S cells contains TF-tagged B7H4 protein, detected with an anti-TF antibody.

#### Products and Services

- Monoclonal Antibody
- Bispecific Antibody
- mRNA-LNP
- Hybridoma Sequencing
- CAR-T/NK
- Stable Cell Line